


## مشخصات محصول

- مناسب برای حجم کاری بالا
- دارای ۱۲ برنامه کاری گوناگون (با زمان و ابعاد دوخت مختلف)
- کیفیت بالای دوخت
- زمان بسیار کوتاه دوخت (بازه زمانی ۲-۴ ثانیه برای هر دوخت)
- قابل استفاده برای کیسه های مختلف خون (از جنس PVC)
- مناسب برای لوله های تا قطر خارجی ۶ میلیمتر
- امکان جداسازی آسان لوله پس از دوخت
- قابلیت اطمینان از کیفیت دوخت و تکرار مداوم آن
- دارای نمایشگر جهت مانیتورینگ لحظه ای
- قابلیت اتصال هند بصورت مجزا به نمونه های رومیزی
- قابلیت اضافه نمودن تعداد ایستگاه کاری متصل با فاصله استاندارد ۷ سانتیمتر تا ۵ ایستگاه
- قابلیت حمل آسان



## Technical Data

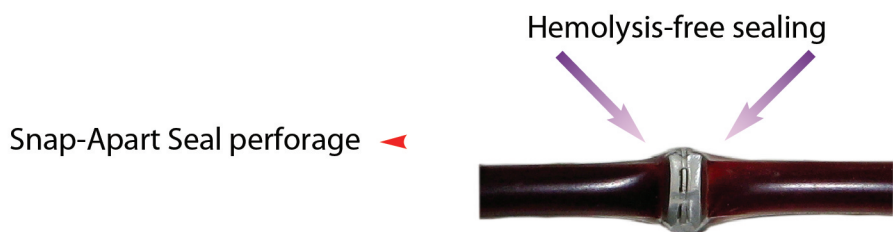

Model	MHSTA	MHSMA	MSHA Portable
Frequency	40.680 MHz	40.680 MHz	40.680 MHz
Sealing Time	1-3 sec (Depended on diameter of tube)	1-3 sec (Depended on diameter of tube)	1-3 sec (Depended on diameter of tube)
RF Output	80 W	80 W	45 W
Electricity Supply	220 V - 50 Hz	220 V - 50 Hz	100-240 V - 47-63 Hz
Dimentions (mm)	70*220*335	70*220*335 each -up to 7 station	80*220*30
Weight	4 kg	4 kg each station	دستگاه: 0.8 kg دسته: 0.3 kg
Sealing Counter	Yes	No	No
Cord Type	Up to 6 mm PVC	Up to 6 mm PVC	Up to 6 mm PVC

### Sealing hand

Sealing hand is available in two models



MHA



MHB

## Save time with a flexible, mobile, handheld tuning sealer

- The MSHA, combines reliability, flexibility and portability to make dependable seals simple and efficient.
- The sealer can be used while the donor is still connected to a disposable tubing set. This versatile device seals polyvinyl chloride (PVC) and ethyl vinyl acetate (EVA) tubing and automatically adjusts to accommodate a range of tubing diameters.
- With a fully charged battery a sealer can make between 1100 and 1500 seals.
- Makes seals with a notch at the center for easy segmenting of tubing.
- Seals wide range of blood bag tubing - good for tubing diameter 2-6 mm.
- Head protector is detachable for sealing head cleaning.
- Canvas carrying bag is available. The bag carries all parts.
- Battery charger is built-in in the unit.
- Battery can charge during sealing.

### RF jack

- Handheld connects to device from this port.

### Power key

- When the power on, lamp turns on


- Light display the status of the sealing process
- The lamp is turned on by pulling the trigger level, and is turned off automatically when a seal is made at tubing.

### Battery charge level & status lamps

- The green LED lamp indicate full battery charge level.
- The yellow LED lamp is turned on when the charger unit plug is plugged into the charging jack.
- The red LED lamp indicate device needs to charge .

### Charging jack

- It is possible to use sealer in both direct power or battery
- Battery can charge during using device

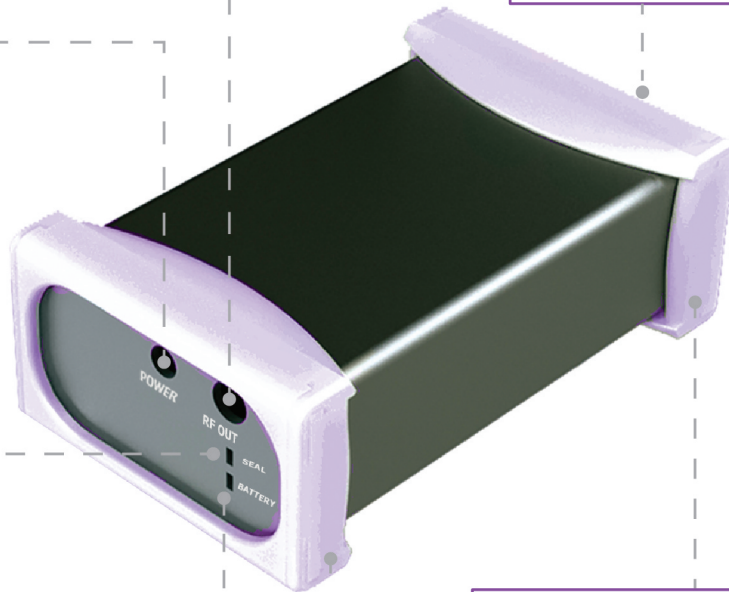

### Silicone protector

- Protection against slip at the time of deployment

## Handheld is available in two modle

