

Teif Gostar Faraz

HID Detector

The HID-2552 is more sensitive detector for analyzing permanent gas and all other gas that TCD detector can analysis but more sensitive

It is better to use combination of both HID and TCD together because high sensitivity of HID get saturated at high concentration but TCD not saturated.



HID-2552

Performance

High speed, High resolution, High reproducibility, High quality,

High sensitive HID

Four filament low noise TCD detector.


Column outlet can be connected to TCD inlet and TCD outlet directly conduct to HID inlet make the instrument best choice for analyzing gaseous component from PPM level to 100%

Application

Permanent gas analyses

All gas analyses

Scientific research

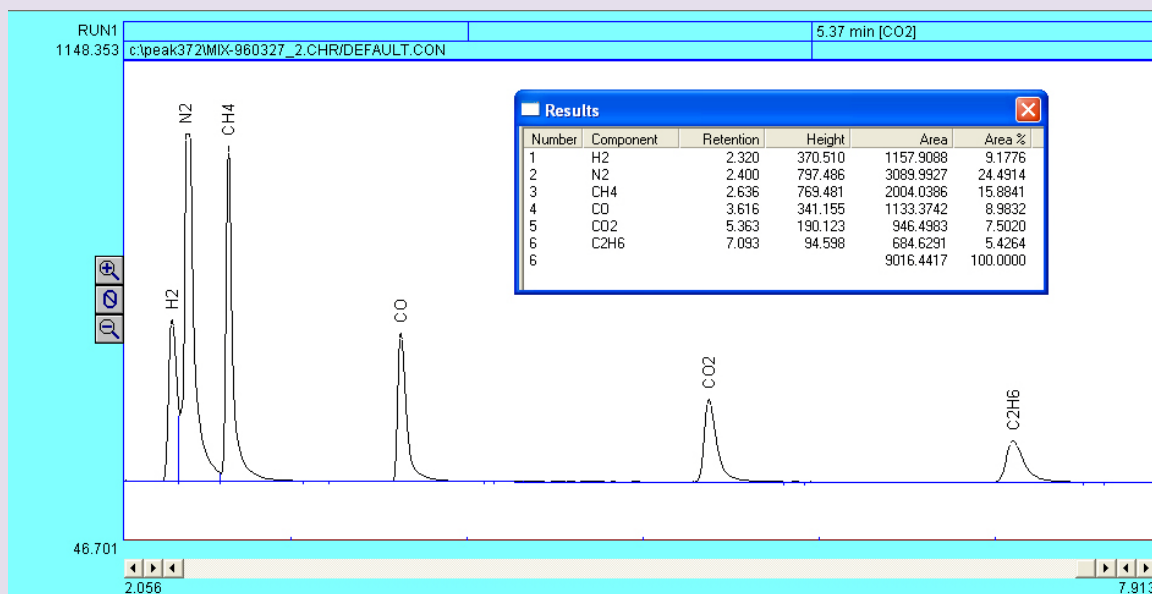

Unit 4, No. 183,
MahtabBldg., TolidGaran Blvd.,
KM 19 Jade Ghadim Karaj,
Tehran 3114859018, Iran

Tel: 0098 21 44336672
Fax: 0098 21 46825835 E
mail: info@irangc.com
Telegram:
[@iranqc95TeifGostarfaraz](https://t.me/iranqc95TeifGostarfaraz)



طیف گستر فراز
طراحی و ساخت گاز کروماتوگراف
تماس: ۴۴۳۳۶۶۷۲